


UnderArm Crutch

User Manual

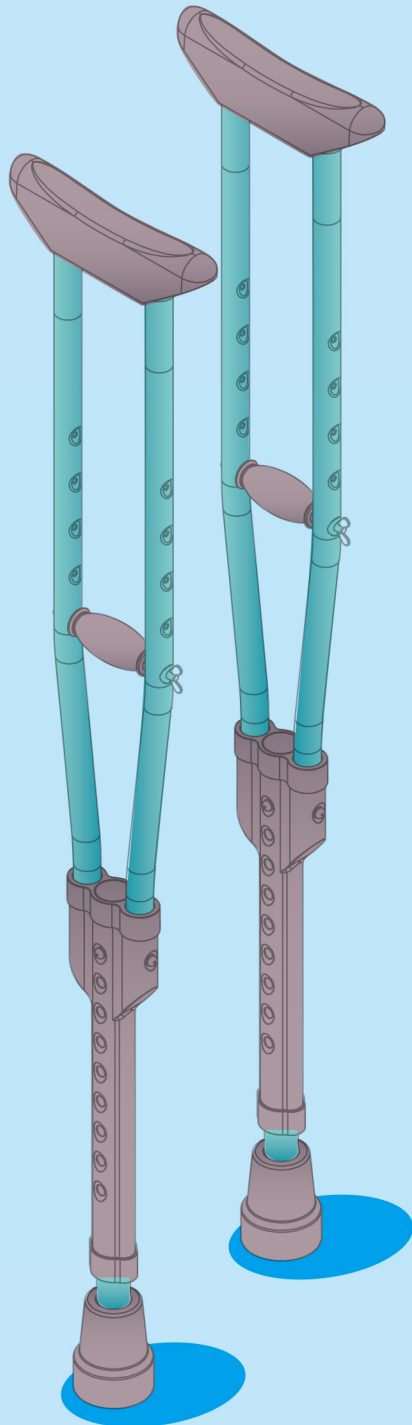


Assembly instructions


•Parts

Main Frame x2 Shaft x 2


Side Frame x2 Tip x 4 Card x 1

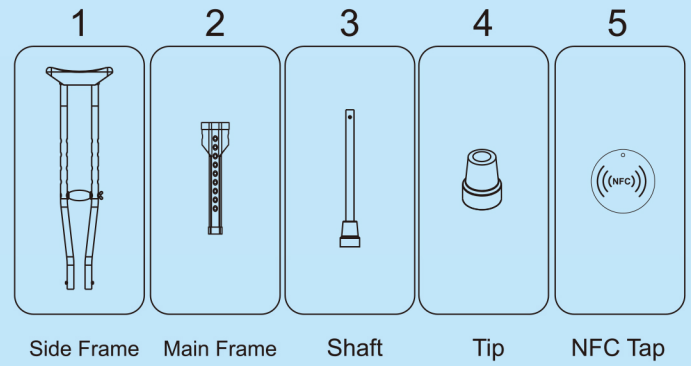

Step 1:

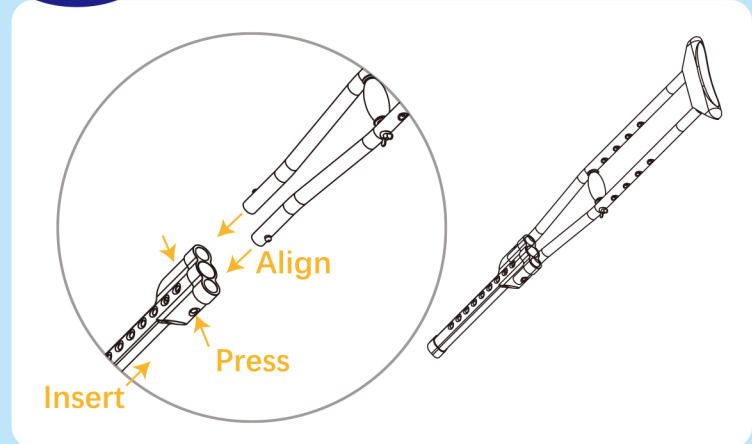


Step 2:

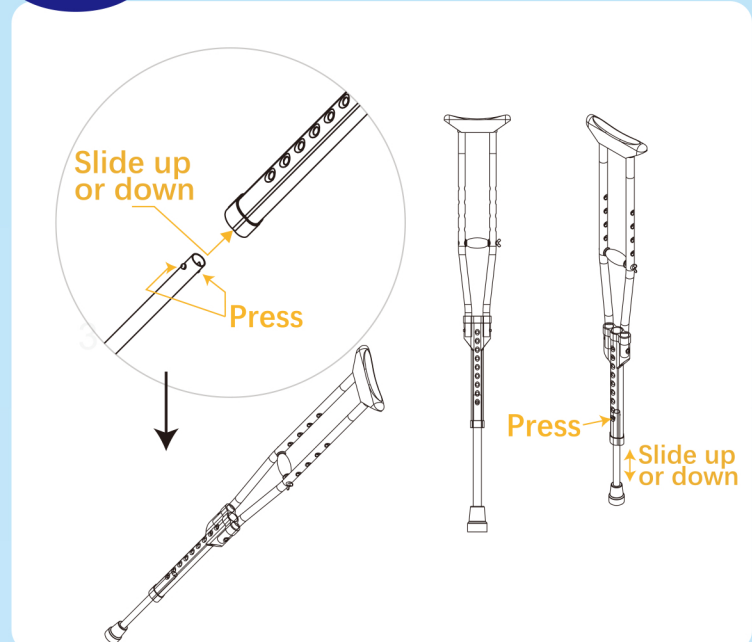

MODEL : HY7528L


Manufacturer: DOUBLE CARE MEDICAL TECHNOLOGY CO., LTD
Address: BUILDING A, SECOND STREET, EAST OF THE SOUTH
SECTION (LIBO WATER SECTION) OF YONGSHI AVENUE, SHIWAN
TOWN, BOLUO COUNTY, HUIZHOU CITY, GUANGDONG

Step 3:

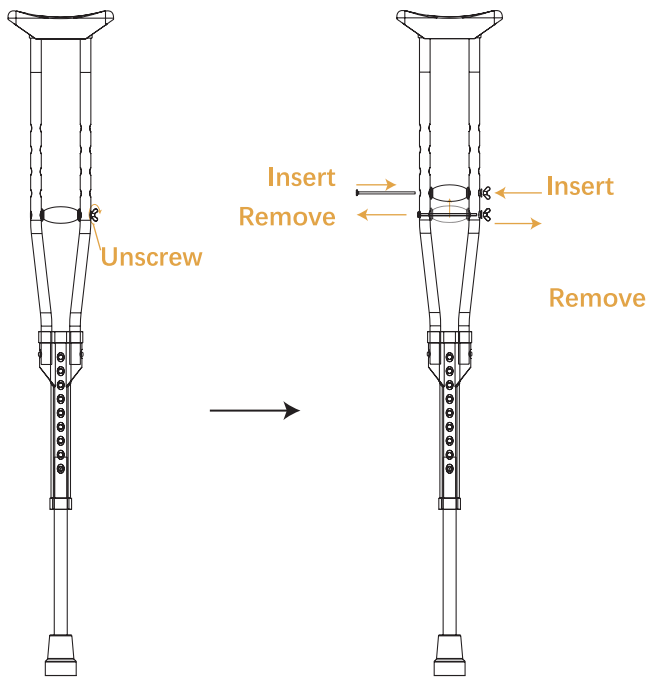

Step 3:

You're All Set!
Perform a final safety check to ensure all connections are secure before using your crutches.

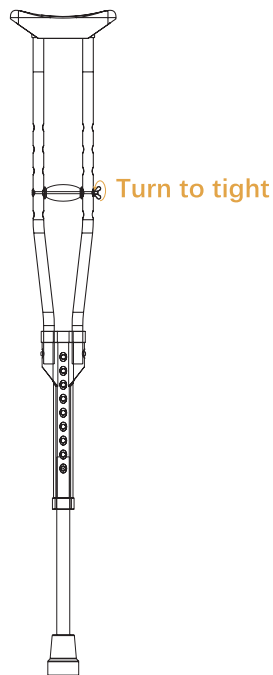

Warranty


This device has been manufactured to the highest standards and carefully inspected before shipping. Enjoy peace of mind with our 1-year warranty and outstanding customer support.; however, it does not cover damage resulting from misuse, neglect, or accidents, nor does it cover normal wear and tear. Non-durable components such as casters, rubber accessories, and grips are also excluded, as they are subject to expected wear and may require periodic replacement. If you have any questions about your DeepTireh device or this warranty, please feel free to contact us.

WARNING


Adjust to correct height. Wrong fit causes pain or falls.

- Never stand on crutches. They are not steps.
- Check for damage before use. If broken—DO NOT USE. Contact seller.
- Get training from a doctor or therapist first.
- No water use (pools, showers). Rust damages parts → injury risk.
- Avoid slippery, soft, or uneven ground. Stay under weight limit.
- Move both crutches together. Don't push one side at a time.
- Lock all buttons fully. Must click and show through holes.
- Inspect rubber tips daily. Replace if worn, torn, or loose.
- Only use approved accessories.
- Do not hang full weight on handgrips.
- Weak or limited users need help/supervision.

DeepTireh Medical Products assumes no responsibility for any damage or injury caused by improper installation, assembly or use of this product.

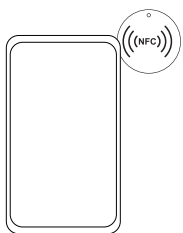

Hold your phone close to this accessory

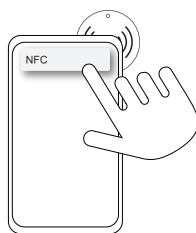

Tap the pop-up link to watch the setup video or view warranty info.

Contact product support , we're always happy to assist!
deeptireh@gmail.com